



**PERCUSSIONAIRE®**  
CORPORATION

# IPV®-2C



**The IPV®-2C** enhances the clinical versatility of IPV® therapy by providing added control, including demand CPAP to reduce work-of-breathing (WOB) while transitioning to Flow Ventilation®, as well as independent CPAP for weaning.

**IPV®-2C**  
K00002-C-OP D

Designed specifically for hospital use and with the versatility to treat most patients, including those in the NICU, the IPV®-2C provides Lung Conditioning Therapy, either invasive or noninvasive, for patients in all hospital units. Flow Ventilation® enables improved clearance of secretions while supporting lung recruitment and gas exchange. The Phasitron® provides a lung protective strategy and also makes the IPV®-2C compatible with conventional ventilators, preventing the need to switch between different devices. It also provides an added advantage of demand CPAP to reduce WOB while transitioning to Flow Ventilation®.

## Accessories

Single Patient use only.

**Phasitron Kit**  
A50095-D



**Phasitron Kit**  
A50095-F



## Technical Specifications




Accessories	Part No.
Phasitron kit	A50095-D or A50095-F

Settings	Specifications
Aerosol flow	Liquid consumption of 1cc per minute
Pulse/interval ratio	Adjustable from 1:1 to 1:5

Operating Conditions	Specifications
Operating range	Temp., 0° C to 49° C (32° F to 120° F), humidity 5%-95%, Elevation 0-9842 ft
Storage and transport range	Temp., -40° C to 60° C (-40° F to 140° F), humidity < 93% non-condensing
Battery type	Multimeter uses ( 2) CR123A battery, Percussionaire part # B13350
Run time	Non-Continuous

Display/Output	Specifications
Pulse Amplitude	Digital display, 0 to 150 cmH2O, accurate to +/- 1cmH2O
Pulse frequency	60-400 pulses per minute
Mean airway pressure	Digital display
Amplitude Bar Graph	Digital display

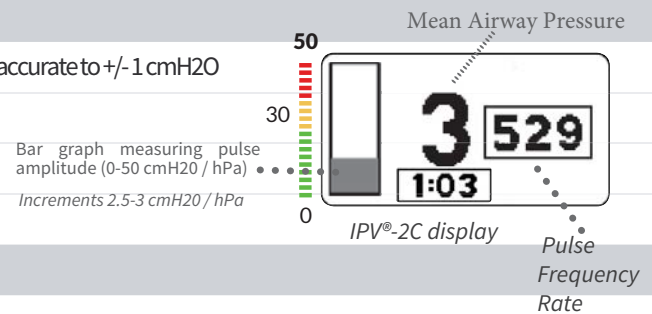

Gas Source	Specifications
Wall gas	50-60 PSI
Flow	25-30 LPM

Filters	Description
Red line filter	Hydrophobic

Dimensions/Weights	Specifications
Dimensions (WxHxD)	6.7" x 9.5" x 9.5"
Weight	4.4 lbs

Service Life	Overhaul every 3 years from date of delivery.
--------------	---



**PERCUSSIONAIRE<sup>®</sup>**  
CORPORATION

130 McGhee Road, Suite 109  
Sandpoint, ID 83864  
800-850-7205 | 208-263-2549  
www.percussionaire.com