




PERCUSSIONAIRE®
CORPORATION

Travel Air®
IMPULSATOR®



The Travel Air® Impulsator® allows patients to benefit from lung conditioning therapy at home. Proven safe and effective, the new Travel Air delivers the same Flow Ventilation, using the Phasitron, as the Percussionaire hospital devices. This self-contained device is easy to use, and is suitable for adult and pediatric patients.

Travel Air®
K00055-US

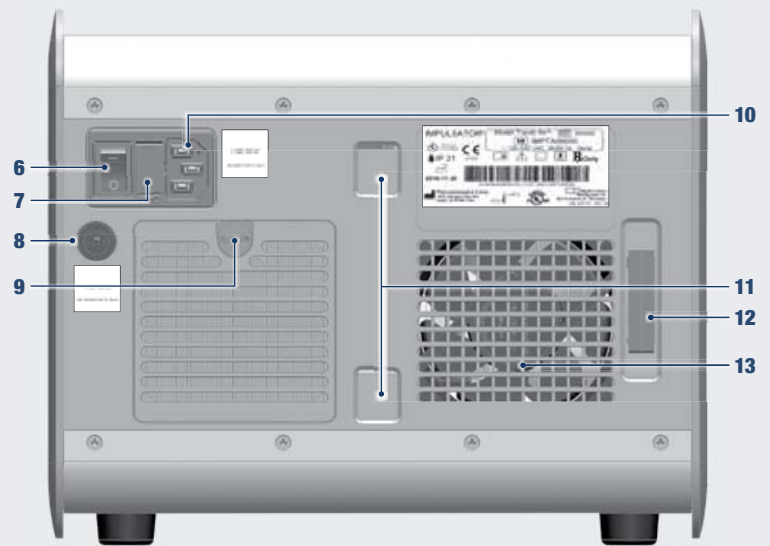

- 1 Multimeter display
- 2 Tubing sockets
- 3 Power indicating LED lights
- 4 Frequency/percussion knob
- 5 Flow rate knob

- 6 Power switch
- 7 Switch fuse
- 8 Thermal fuse
- 9 Release latch for cooling air filter
- 10 Power cord socket

- 11 Handle or stand clamp mounting sockets
- 12 Patient air inlet filter
- 13 Cooling air outlet




Phasitron Kit
A55108

Single Patient use only.

Technical Specifications

Accessories	Part No.
Phasitron kit	A55108
Settings	Specifications
Aerosol flow	Liquid consumption of 1cc per minute
Percussion frequency	up to 600 cycles/minute
Pulse/interval ratio	1:2 - 1:5 Automatic
Operating Conditions	Specifications
Operating range	Temp., 5° C to 32° C (41° F to 90° F), humidity 5%-95%, Elevation 0-9842 ft
Storage and transport range	Temp., 10° C to 60° C (50° F to 140° F), humidity <93% non-condensing
Power supply	100-240 VAC, 50/60Hz, 3 amp
Run time	30 minutes ⚠ Caution Do not exceed 30 minutes, allow 10 Minute cool down before restart.
Display/Output	Specifications
Airway pressure	Digital display, 0 to 150 cmH2O, accurate to +/- 1cmH2O
Pulse frequency	Digital display, 50-999 cycles/minute.
Usage Timer	Digital display
Total operating time duration	Digital display
Condition	Indication
Airway pressure over 90-100 cmH2O	Audible indicator
Compressor overload	Audible indicator/shutting down of the Travel Air®
Compressor over temperature	Shutting down of the Travel Air®
Compressor	Specifications
Pressure	40 PSI
Flow	25 LPM
Filters	Description
Patient air filter	Compressor entrance filter felt
Red line filter	Monitoring line filter
Cooling air filter	80mm filter intake
Dimensions/Weights	Specifications
Dimensions (WxHxD)	10.5" x 10" x 10"
Weight	19.8 lbs
Service life	6 years from date of delivery

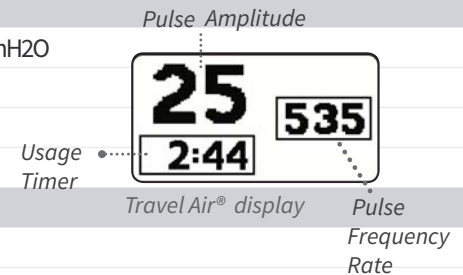

PERCUSSIONAIRE®
CORPORATION

130 McGhee Road, Suite 109
Sandpoint, ID 83864
800-850-7205 | 208-263-2549
www.percussionaire.com